




Akja Windshield Cover



Fixation at Akja

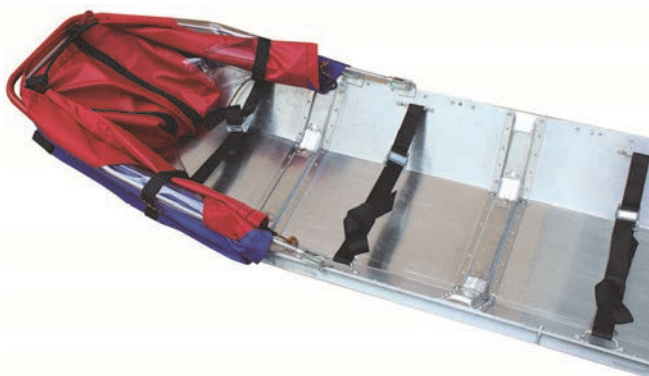

Raise and tighten the cover



Provided parts „Windshield cover“

Roomy, foldable wind protection cover that provides effective protection for the injured in the Akja against airstream and whirled-up snow - especially, if the Akja is towed behind a snowmobile.

Due to the large window areas, the injured feels unobstructed.

Easy loading of the injured.

Due to the flexible fixing brackets the windshield cover could be completely folded and could always remain at the Akja, even if the windshield cover is not in use. In this case, a vacuum mattress or injured rescue bag will hide the cover.

The cover is **attached to the Akja** with the provided quick fastener fork head system and straps.

Completely **manufactured** from high-strength, water resistant nylon-fabrics. This **high-grade material** ensures longevity, maximum wear-resistance, pleasant handling and does not become stiff even under very cold temperatures.

The transparent **window material** used in the field of offshore yachting withstands heavy duty conditions.

Suitable for **all TYROMONT Akja** models.

Art. No.	93063	
Dimensions	Length:	1.200 mm
	Width:	580 mm
	Height:	700 mm
Weight	2,0 kg	
Comments	Ideal solution together with Akja Coupling Device. (Art. 93061)	